

CORDLESS ENDO MOTOR

CHANGZHOU BOMEDENT MEDICAL TECHNOLOGY CO., LTD



Add: N9. Changyang Road West Taihu Science&Technology Industry
Park Changzhou, Jiangsu, China



E-mail: anita.wu@bomedent.com



Tel: 0086-0519-89881606

WISMY

www.bomedent.com



Changzhou Bomedent Medical Technology Co., Ltd.

Changzhou Bomedent Medical Technology Co., Ltd. is a high-tech enterprise integrating R & D, design, production, maintenance and sales. The main products are root canal therapy instrument, root canal length measuring instrument and its peripheral products. There are many clinics, and most of the high-end dental medical equipment is imported. The domestic dental medical equipment does not match the growing demand for oral medical treatment in China. Dental treatment is generally more expensive. In order to change the status quo, reduce the cost of dental treatment, and improve the level of dental treatment of Chinese people, our company will closely integrate the group's purpose: establish a new domestic dental medical ecological environment, research and development, production, sales of high-end dental medical equipment, the company has a A research and development team with more than ten years of experience in dental product research and development. Our team will uphold the concept of technology as the core and quality as life, and continue to develop and innovate. We will provide you with cost-effective products and meticulous after-sales service.

CORDLESS ENDO MOTOR WITH INTEGRATED APEX LOCATOR

WISMY

TYPE	PORTABLE/CORDLESS
MEMORY	10 MEMORY SETTING
ROTATION	CLOCKWISE/ANTI-CLOCKWISE/ RECIPROCATING
RECOPROCATING	ROTATION ANGLE ADJUSTABLE
APEED/TORQUE	800 RPM MAX/4.0N.CM MAX
BATTERY	RECHARGABLE/REPLACEABLE
CHARGING	MICRO USB
DISPLAY	OLED
AVAILABLE CONTRA ANGLE	16:1
APEX LOCATOR CONTROL MOTOR WORKING MODE/TORQUE CONTROL	AUTO START IN THE ROOT CANAL AUTO REVERSE WHEN REACHING APICAL FORAMEN AUTO REVERSE WHEN EXCEEDING THE TORQUE LIMIT VALUE

FOUR WORKING MODES

1

APEX
LOCATOR
CONTROL
MOTOR

2


MOTOR
APEX
INDEPENDENT

3

ONLY
APEX
LOCATOR

4


ONLY
MOTOR

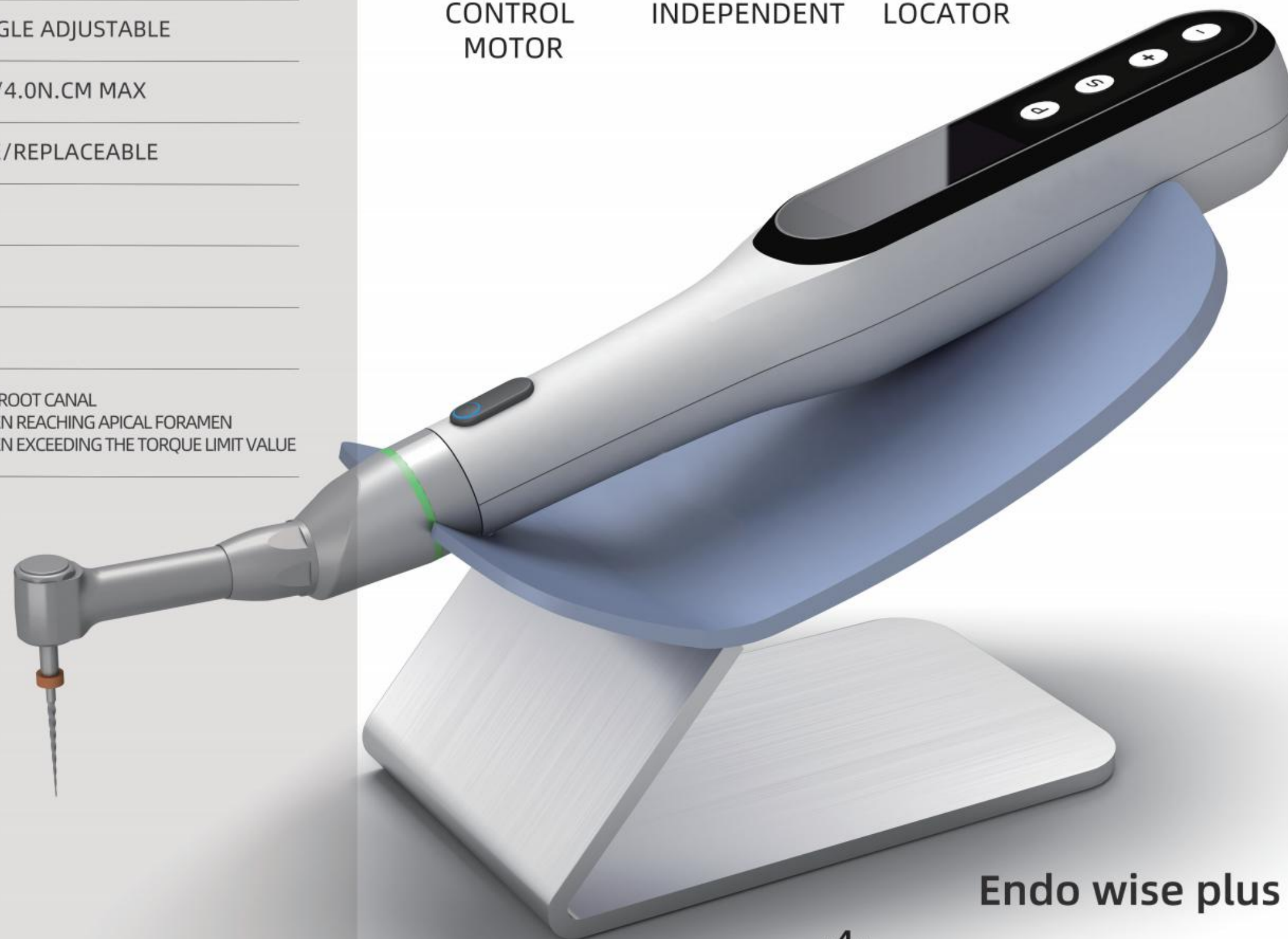


ENDO MOTOR
+
APEX LOCATOR

Endo wise plus

Apex locator control motor working mode:
Auto start/Auto stop/Auto reverse

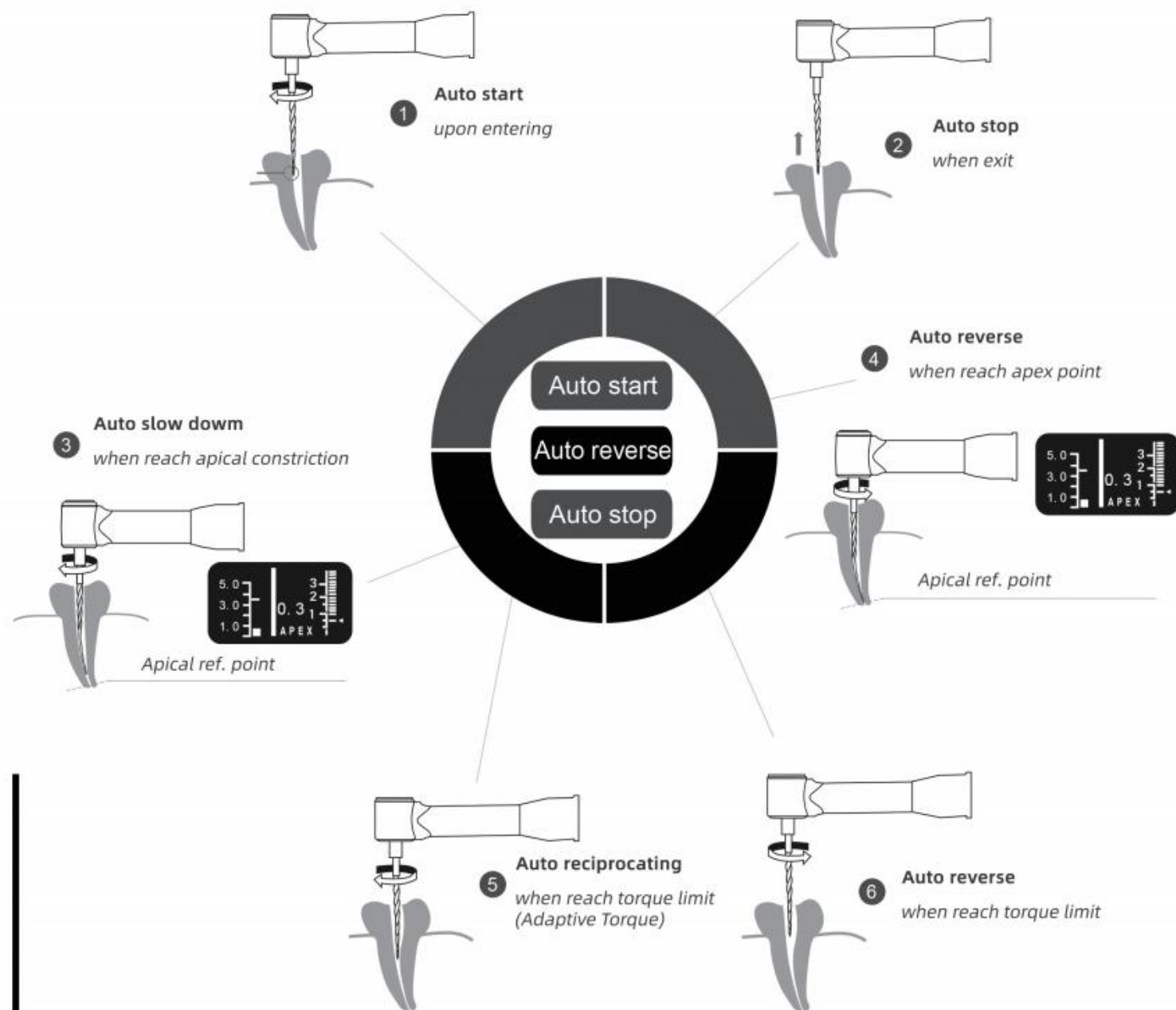

ENDO MOTOR WITH INTEGRATED APEX LOCATOR

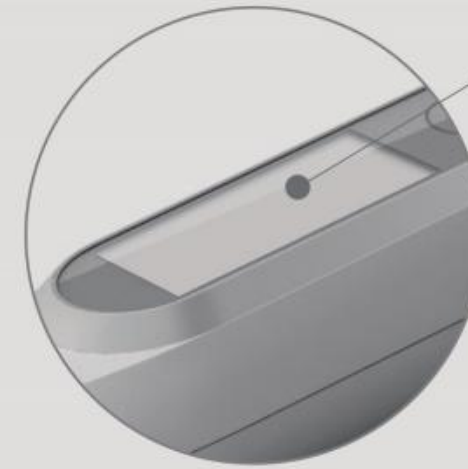

- Large OLED display
- Real-time Torque & WL readings
- Adjustable WL reference position



1. Auto start in the root canal, Auto reverse at apex, Auto stop
2. Reciprocating rotation angle adjustable
3. Excellent power control circuit, long working time

CORDLESS ENDO MOTOR




Endo wise


FEATURES

- 1. Multiple reciprocating modes
- 2. Large OLED Screen, Long working time
- 3. Accurate torque control

TYPE	PORTABLE/CORDLESS
MEMORY	10 MEMORY SETTING
ROTATION	CLOCKWISE/ANTI-CLOCKWISE/RECIPROCATING
TORQUE CONTROL	AUTO REVERSE/AUTO STOP
BATTERY	RECHARGABLE/REPLACEABLE
CHARGING	MICRO USB
DISPLAY	OLED
AVAILABLE CONTRA ANGLE	16:1

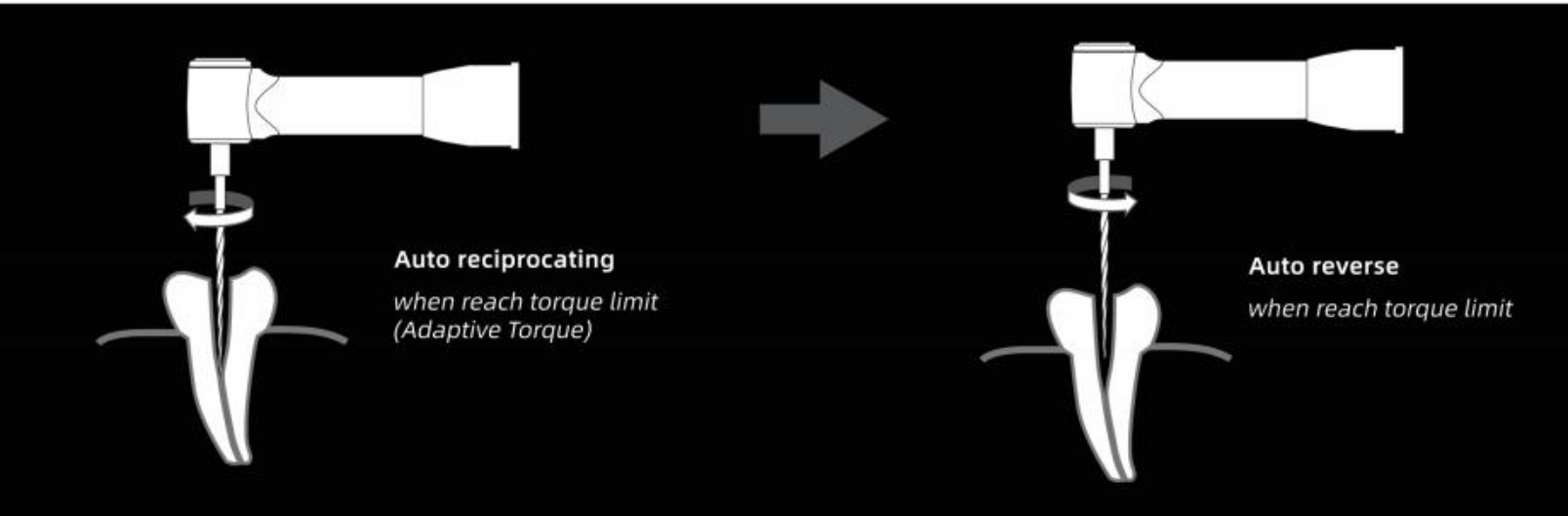

Push button, convenient

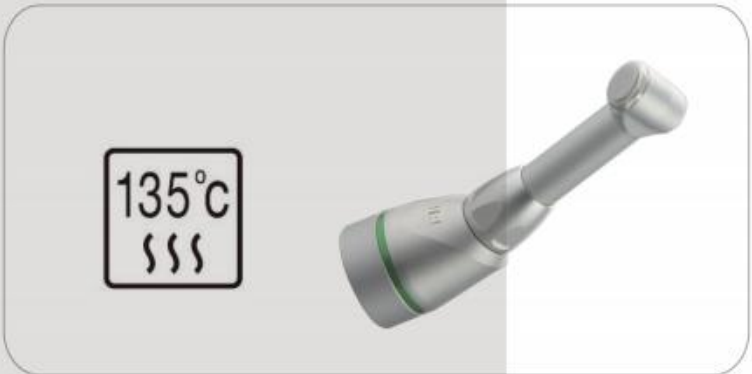

Can be sterilized by high temperature and high pressure



Compatible with most NI-TI file on market , flexible

APEX LOCATOR WISPEX



ACCURATE


SENSITIVE



- High definition TFT screen.
- Multi frequency measurement technology guarantees high accuracy.
- Multi color regions.
- Segmented display of the file in root canal.
- Working length adjustable.



Not affected by tooth types

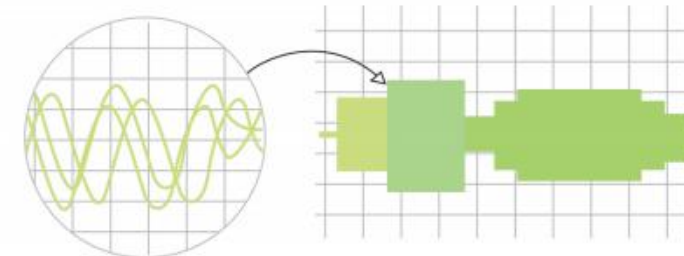

Blood or pulp

Multi frequency measurement technology.
Not affected by the blood or pulp in the root canal.